

# Reciprocal Space (a Gentle Introduction)

Paul D. Boyle  
Department of Chemistry  
University of Western Ontario  
London, ON

# Introduction

- Notions of **reciprocal space** and the **reciprocal lattice** are indispensable for understanding X-ray Crystallography
- Diffraction by a lattice gives a lattice as a diffraction pattern
- There is a precise mathematical relationship between the original diffracting lattice and the resulting pattern
- However, this is a (mostly) non-mathematical introduction to help you qualitatively understand the concepts of reciprocal space and the reciprocal lattice

# Direct Space & Reciprocal Space

- We live in direct space
- Distances and orientations between isolated objects
- Reciprocal space is a “spatial frequency” space (e.g. number of Tim Horton’s per kilometre)
- In NMR time and frequency are related by a Fourier transform (units: time  $t$  and frequency  $t^{-1}$ )
- In X-ray Crystallography direct space and reciprocal space are related by Fourier transform (units: distance  $\text{\AA}$  and spatial frequency  $\text{\AA}^{-1}$ )
- For every direct space lattice there is a corresponding reciprocal space lattice

# Bragg's Law

- $n\lambda = 2d\sin\theta$
- Rewrite as:
- $\sin\theta = 0.5n\lambda(1/d)$
- Reciprocal relationship between diffraction angle,  $\theta$ , and the  $d$  spacing
- The smaller the  $d$  spacing, the higher the diffraction angle

# Reciprocal Quantities


- $1/d = d^*$
- $d^*$  is a reciprocal quantity and typically has units of  $\text{\AA}^{-1}$
- The '\*' in general means a reciprocal space quantity
- A direct space unit cell with parameters  $a, b, c, \alpha, \beta, \gamma$  has a corresponding reciprocal unit cell:  $a^*, b^*, c^*, \alpha^*, \beta^*, \gamma^*$
- A reciprocal lattice vector:
  - Has the form  $\mathbf{d}^* = h\mathbf{a}^* + k\mathbf{b}^* + l\mathbf{c}^*$
  - The values  $h\ k\ l$  are the components of a reciprocal lattice vector
  - The vector  $\mathbf{d}^*$  is perpendicular to the set of  $(hkl)$  lattice planes
- There are exact mathematical relationships which relate the direct space and reciprocal space unit cell parameters
- We can calculate the **reciprocal metric tensor,  $\mathbf{G}^*$** 
  - $\mathbf{G}^* = \mathbf{G}^{-1}$  and  $\mathbf{G} = \mathbf{G}^{*-1}$
  - $d^{*2} = \overline{\mathbf{h}}\mathbf{G}^*\mathbf{h}$  where  $\mathbf{h} = [h\ k\ l]$

# Reciprocal Relationships

- The relationships between the direct axes and reciprocal axes is strictly reciprocal
- Any statement about the two lattices remains true if you simply replace all starred (\*) quantities by unstarred quantities and vice-versa
- $a^* \perp bc(\text{face})$  and  $a \perp b^*c^*(\text{face})$
- Any direct axis has as family of reciprocal lattice planes which are perpendicular to that axis
- Conversely, any reciprocal axis has a family of direct lattice planes which are perpendicular to that axis

# X-ray Diffraction Patterns

- The X-ray diffraction pattern is the reciprocal lattice of a crystal's direct lattice
- Referred to as the **intensity weighted reciprocal lattice**
- Diffraction maxima are reciprocal lattice points
- Intensity distribution of diffraction pattern is related to the electron density distribution in the crystal

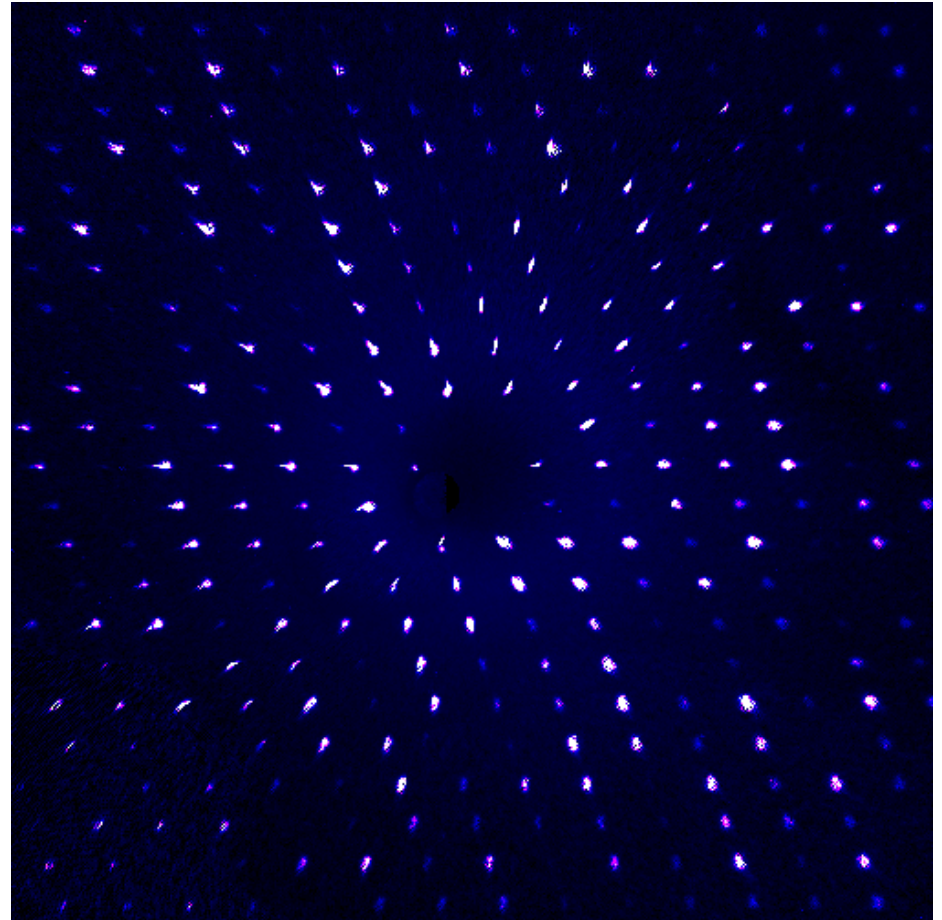


# Reciprocal Lattice Points

- Are designated by their Miller index,  $hkl$
- Assigning  $hkl$  values to the reciprocal lattice points is called ***indexing the crystal*** or ***indexing the diffraction pattern***
- Reciprocal lattice points represent the diffraction from a **set of planes** designated by the  $hkl$  value and have a corresponding  $d^*$  value
- Normal to the set of planes and therefore represent a direction in reciprocal space




# Graphical Construction of Reciprocal Lattice from Direct Space Lattice

- For a set of planes in direct space, we draw a vector normal to these planes
- Terminate the vector at a distance  $1/d$
- For a given lattice row:
- $d^*(nh,nk,nl) = nd^*(hkl)$
- Graphic:  
[http://www.xtal.iqfr.csic.es/Cristalografia/parte\\_04-en.html](http://www.xtal.iqfr.csic.es/Cristalografia/parte_04-en.html)

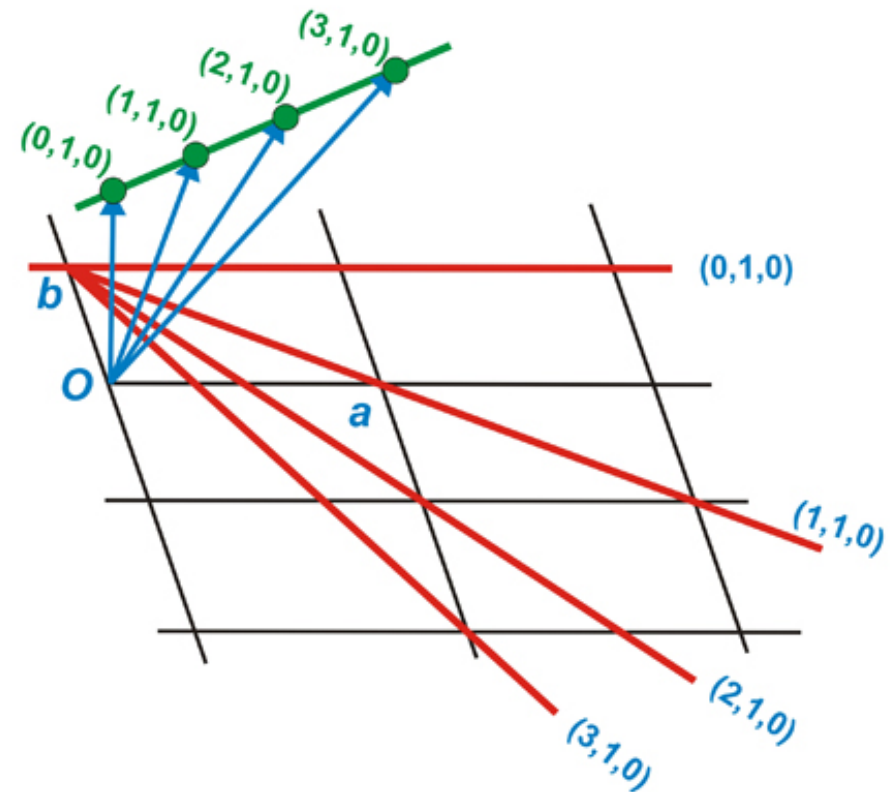


# Indexing a Diffraction Pattern

- Synthesized reciprocal lattice layer (hk0) from an actual crystal
- Vertical axis has closer packed reciprocal lattice points
- Vertical axis has larger direct space unit cell parameter

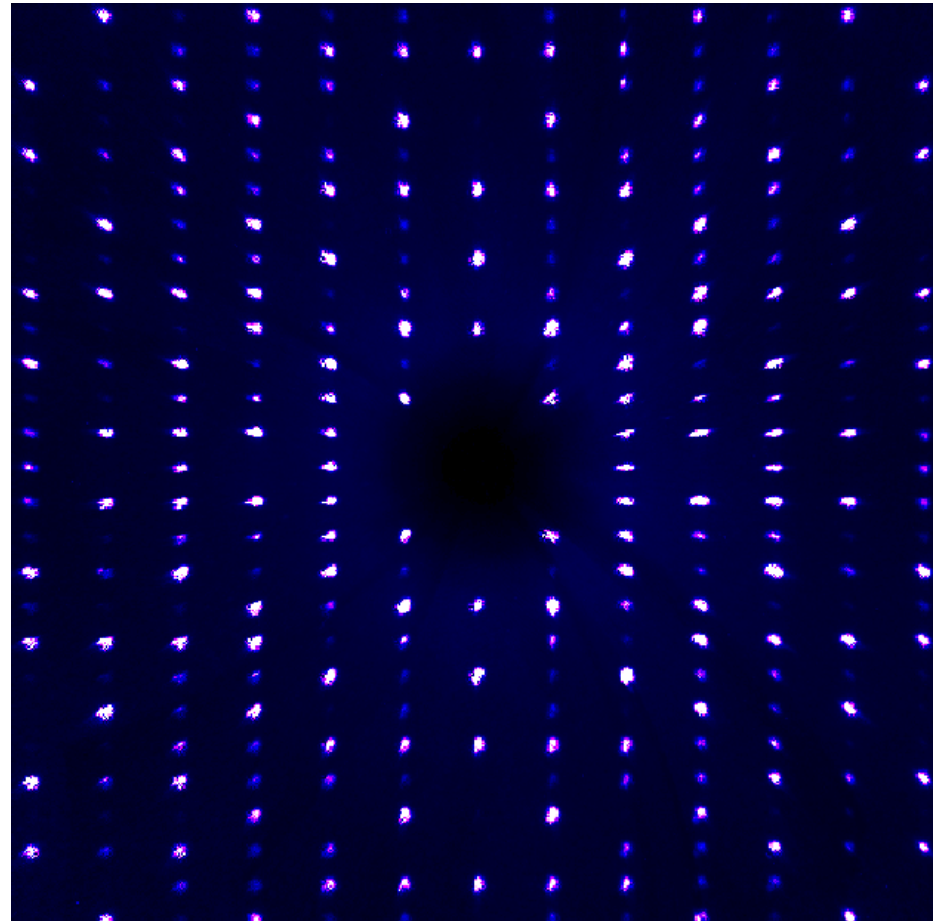


# Indexing a Diffraction Pattern

- First assign the lattice directions
- Notice there are systematic absences along the  $h00$  and  $0k0$  reciprocal axes
- Indicative of two screw axes (translational symmetry elements)

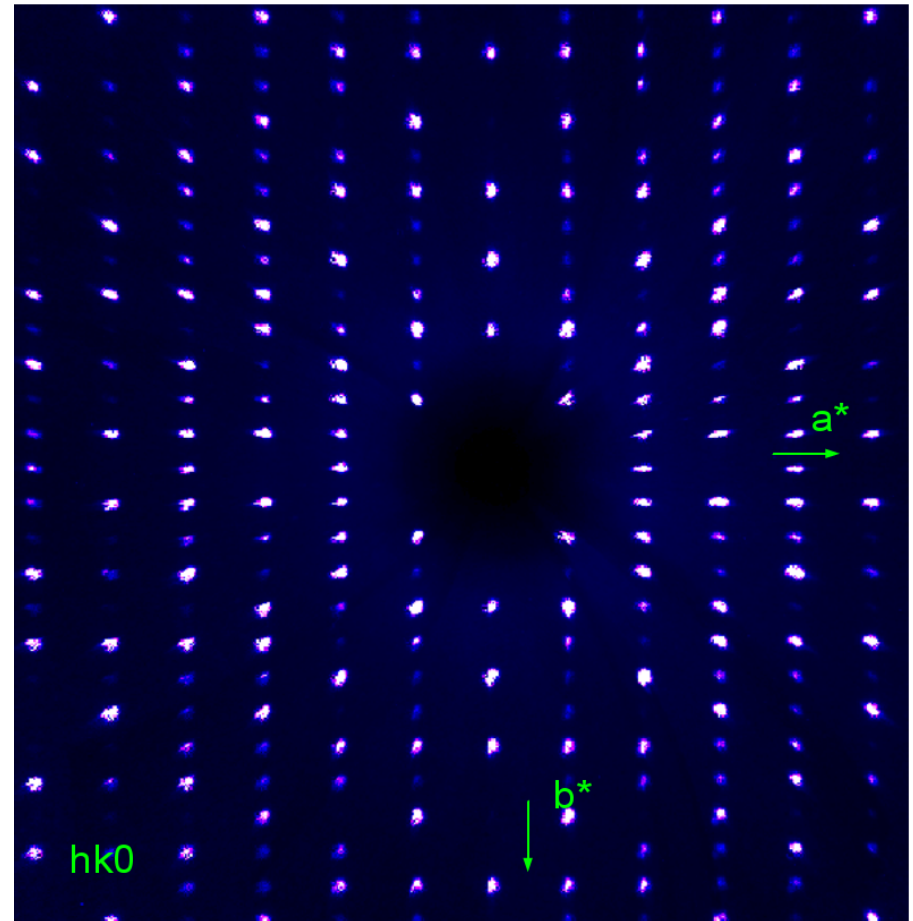


# Indexing a Diffraction Pattern

- Assign  $hkl$  values to each reciprocal lattice point
- Use Bragg's Law to calculate the interplanar spacing associated with each reciprocal lattice point
- Measure angle between  $a^*$  and  $b^*$  to obtain  $\gamma^*$
- Repeat process with other zero layers ( $0kl$  and  $h0l$ )

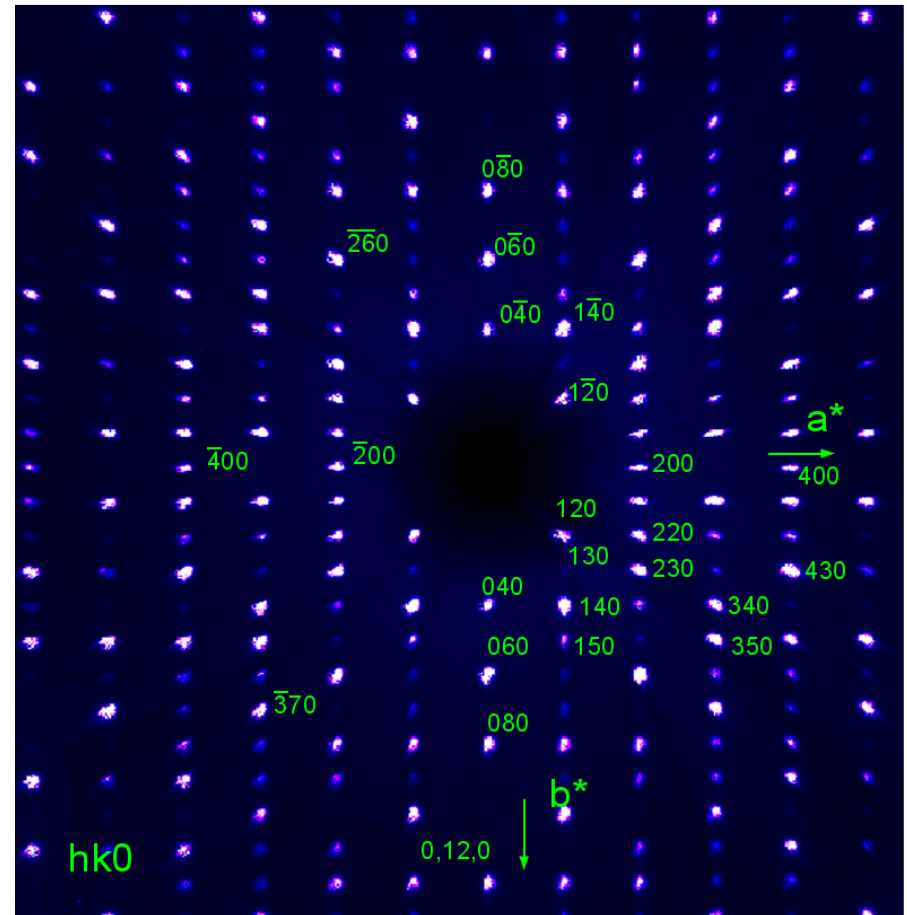

# How to think about this


- Each reciprocal lattice point represents both a direction and d spacing
- With each reciprocal lattice point measured, we are “sampling” the electron density with certain spatial frequency in a given direction

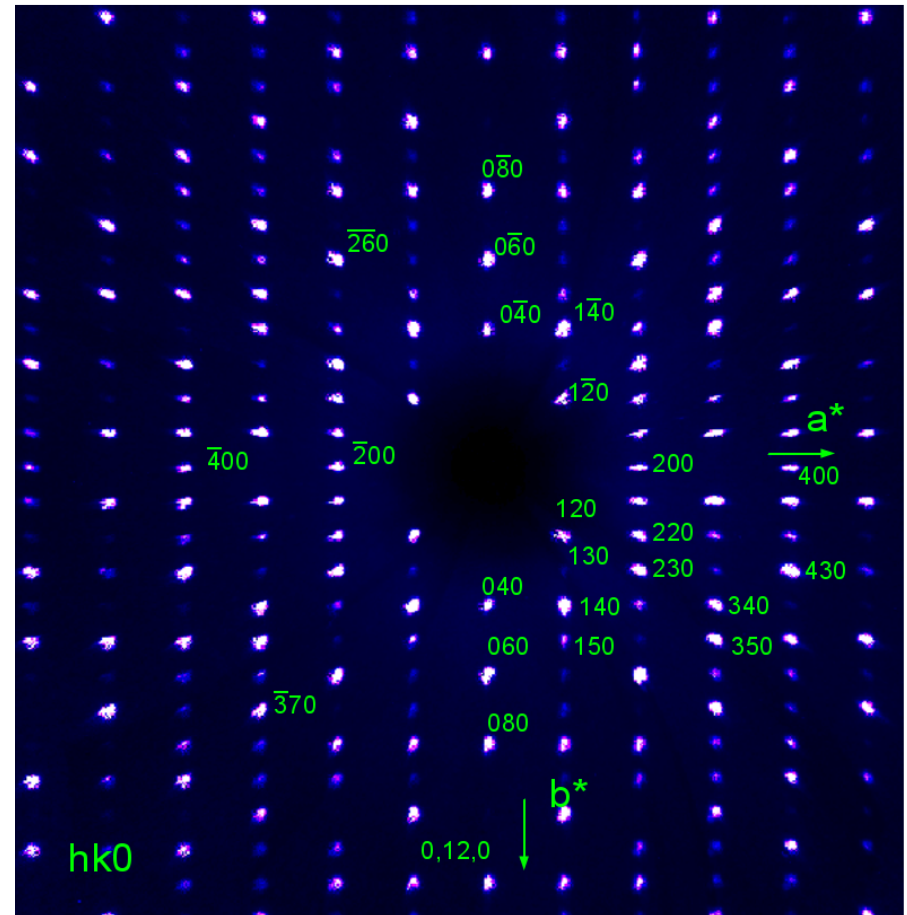




# The Swiss Cheese Analogy

- We want to map where all the holes are in a block of Swiss cheese
- We (virtually) slice the block using various thicknesses and at various orientations
- We then take these slices and use them to map the size and shapes of the all the holes in the block

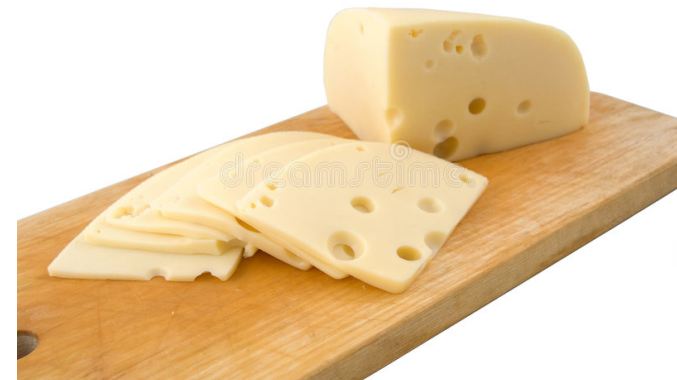


# Resolution of Our Mapping

- Slicing our cheese every 10 mm will cause us to miss some of the smaller holes in the cheese
- We make finer and finer slices to map even the smaller holes within the cheese
- Why not just use all fine slices rather than both low and high resolution slices?
- Analogy breaks down at this point
- In X-ray, we model the electron density as a Fourier series and we need both the low resolution data and high resolution data

Slice Width	Reciprocal Units	Low Resolution
Every 10.0 mm	$1/10$ or $0.1 \text{ mm}^{-1}$	
Every 5.0 mm	$1/5$ or $0.2 \text{ mm}^{-1}$	
Every 2.0 mm	$1/2$ or $0.5 \text{ mm}^{-1}$	
Every 1.0 mm	$1/1$ or $1 \text{ mm}^{-1}$	
Every 0.5 mm	$1/(1/2)$ or $2.0 \text{ mm}^{-1}$	
		High Resolution

# Resolution in Reciprocal Space


- The higher the diffraction angle, the finer the slice we are using to sample our crystal's electron density
- Diffraction condition only allows us to sample the electron density distribution at certain spatial frequencies (Bragg's Law)
- We need to collect both high and low resolution data

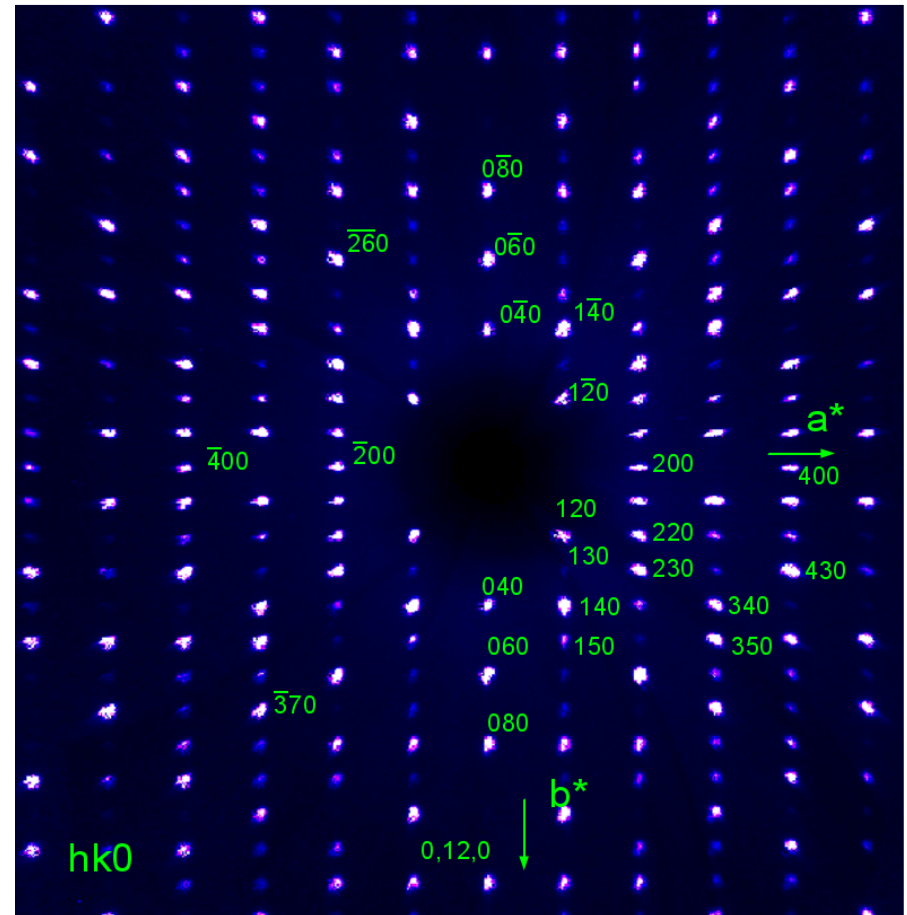




# Ewald Construction

- Graphical depiction of Bragg's Law
- Circle has radius of  $1/\lambda$ , centre at C such that origin of reciprocal lattice, O, lies on circumference
- XO is the X-ray beam, P is the reciprocal lattice point (in this case the 202 reflection)
- OP is the reciprocal lattice vector ( $\mathbf{d}^*$ ) and is normal to the (202) set of planes [aka the Scattering Vector]
- Angle OBP is  $\theta$ , the Bragg angle
- Angle OCP is  $2\theta$
- CP is the direction of the diffracted beam
- BP is parallel to the set of (202) planes
- **Any time a reciprocal lattice point falls on the circumference, Bragg's Law is fulfilled**

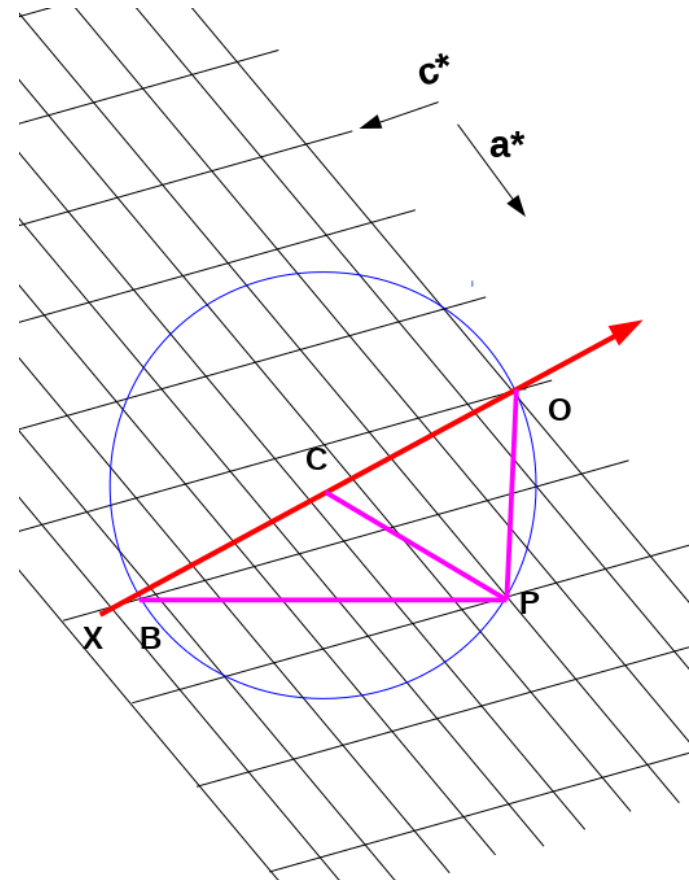

# Ewald Sphere

- 2D Ewald construction can be generalized to 3D to generate the “Ewald Sphere” (also called the “Sphere of Reflection”)
- Anytime a reciprocal lattice point is on the surface of the sphere Bragg’s Law is fulfilled
- Experimentally, we rotate the crystal (lattice) to bring a greater number of reciprocal lattice points pass through the surface of the sphere
- Image shows the detector slicing through part of the Ewald sphere and all the lattice points which were laying on the surface of the sphere

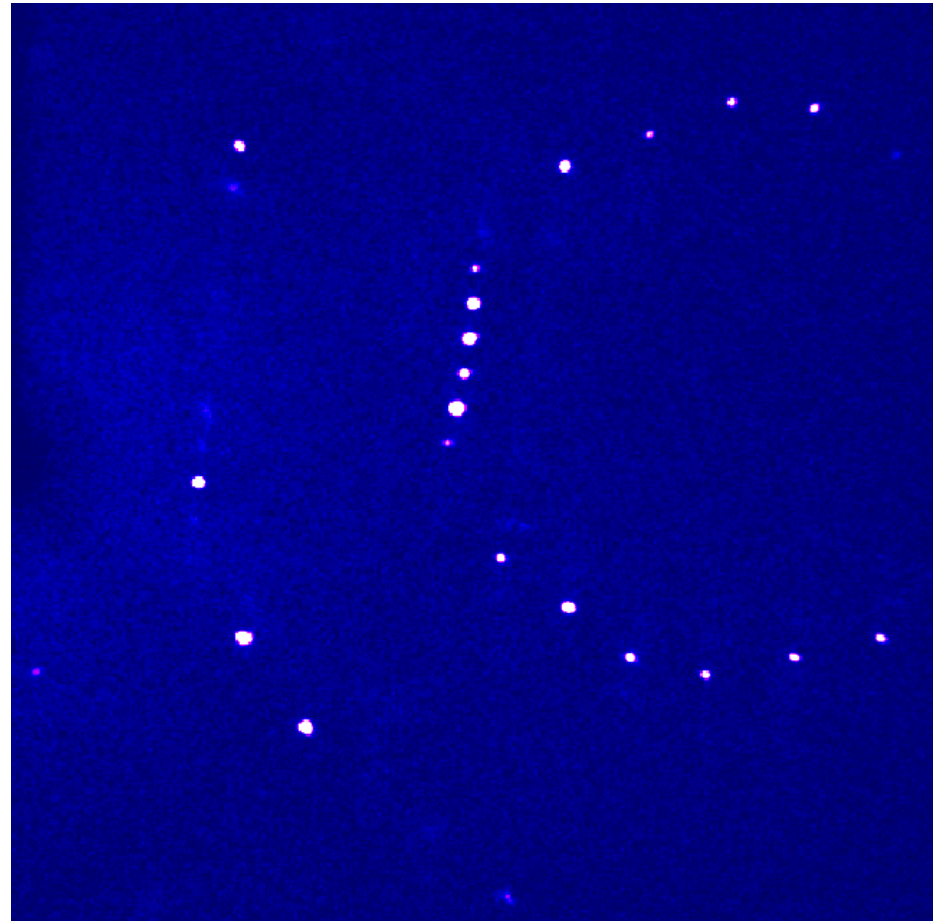

# Ewald Spheres and Limiting Spheres

- Ewald sphere has a diameter of  $2/\lambda$
- Every reciprocal lattice point within that distance can be brought into diffracting position
- Limiting sphere has a radius of  $2/\lambda$
- The total number of reciprocal lattice points within the limiting sphere is approximated by
- $N \approx 33.5(V_{\text{cell}} / \lambda^3)$

# Limiting Spheres of Common Radiations

- $N_{\text{MoK}\alpha} \approx 33.5V_{\text{cell}} / 0.71073^3 = 93.3V_{\text{cell}}$
- $N_{\text{CuK}\alpha} \approx 33.5V_{\text{cell}} / 1.54178^3 = 9.14V_{\text{cell}}$
- Normally, we don't collect all reflections within the limiting sphere. In practice, we pick some maximum value of  $\theta$
- $N_{\theta(\text{max})} \approx (33.5 / \lambda^3)V_{\text{cell}} \sin^3\theta_{\text{max}}$
- You will always get more data with a shorter wavelength

# Wavelength Imposed Limits

- Maximum value of sine function = 1.0
- Imposes certain limits on the X-ray experiment
- Shorter wavelengths allow collection of more data points out to higher resolution

Quantity	CuK $\alpha$	MoK $\alpha$
$\lambda$	1.54178 Å	0.71073 Å
$(\sin\theta/\lambda)_{\max}$	0.648 Å <sup>-1</sup>	1.407 Å <sup>-1</sup>
$d_{\min}$	0.771 Å	0.355 Å
Resolution Limit ( $0.92d_{\min}$ )	0.71 Å	0.33 Å

# Practical Considerations for Data Collection

- Long axes give densely packed reciprocal lattice rows
- Integration is better if peaks aren't overlapping
- Choose minimum crystal to detector distance as:  
$$DX(\text{mm}) = 2 * \text{longest primitive axis } (\text{\AA}) [\text{MoK}\alpha]$$
$$DX(\text{mm}) = 1 * \text{longest primitive axis } (\text{\AA}) [\text{CuK}\alpha]$$
- For non-merohedrally twinned samples, move the detector back even farther

# Determining the Laue Symmetry from a Diffraction Pattern

- The symmetry of a diffraction pattern follows its Laue symmetryWhat does this mean?
- The intensities of (point) symmetry equivalent diffraction maxima are about the same.
- Bruker software: calculates a quantity, “ $R_{\text{sym}}$ ” to give a qualitative indication of Laue symmetry
- An  $R_{\text{sym}}$  value close to zero (  $< 0.07$ ) is a good indicator of symmetry equivalence
- Use point group symmetry operators to see which sets of hkl values should be equivalent

# Equivalent hkl Examples

- **Monoclinic (2/m)**

- $I_{hkl} = I_{h\bar{k}l} = I_{\bar{h}k\bar{l}} = I_{\bar{h}\bar{k}l}$
- $I_{\bar{h}kl} = I_{h\bar{k}\bar{l}} = I_{hkl\bar{l}} = I_{h\bar{k}l\bar{l}}$
- However:  $I_{hkl} \neq I_{\bar{h}kl}$

- **Tetragonal (4/m)**

- Can be used to distinguish 4/m from 4/mmm
- $I_{hkl} = I_{\bar{k}hl} = I_{\bar{h}\bar{k}l} = I_{k\bar{h}l} = I_{hkl\bar{l}} = I_{\bar{k}h\bar{l}} = I_{\bar{h}kl\bar{l}} = I_{kh\bar{l}}$



# Caveat: Point Symmetry in Reciprocal Space!!

- You might assume that the point symmetry operators for reciprocal space are identical to those of direct space
- Nope!
- If  $\mathbf{S}$  is a point group symmetry operator matrix in direct space, the corresponding operator in reciprocal space,  $\mathbf{V}$ , is:
- $\mathbf{V} = \overline{\mathbf{S}}^{-1}$
- $\mathbf{V}$  is the inverse transpose of  $\mathbf{S}$
- This becomes important in trigonal and hexagonal groups

# Example 3-fold axis along [001]

*Direct Space*

$$\begin{pmatrix} 0 & \bar{1} & 0 \\ 1 & \bar{1} & 0 \\ 0 & 0 & 1 \end{pmatrix}$$

*Reciprocal Space*

$$\begin{pmatrix} \bar{1} & \bar{1} & 0 \\ 1 & 0 & 0 \\ 0 & 0 & 1 \end{pmatrix}$$

# Experimental Determination of Space Group

- Space groups are determined primarily through the examination of systematic absences in the diffraction pattern
- Systematic absences arise from the presence of translational symmetry elements
  - Non-primitive lattice centrings
  - Screw axes (rotation with translation)
  - Glide planes (reflection with translation)


# Systematic Absences due to Non-Primitive Lattices

- Non-primitive lattices exhibit systematic absences in the general hkl class of reflections

Centring	Absence Condition for hkl reflections
A	$k+l = \text{odd}$
B	$h+l = \text{odd}$
C	$h+k = \text{odd}$
F	$k+l = \text{odd},$ $h+l = \text{odd},$ $h+k = \text{odd}$
I	$h+k+l = \text{odd}$

# Screw Axis Absences

- Screw axes affect the classes of axial reflections:  $h00$ ,  $0k0$ , and  $00l$
- The type of screw axis is determined by examining the pattern of the absence
- Example: In this figure there is a  $2_1$  axis parallel to  $b^*$
- $0k0$ :  $k = \text{odd}$

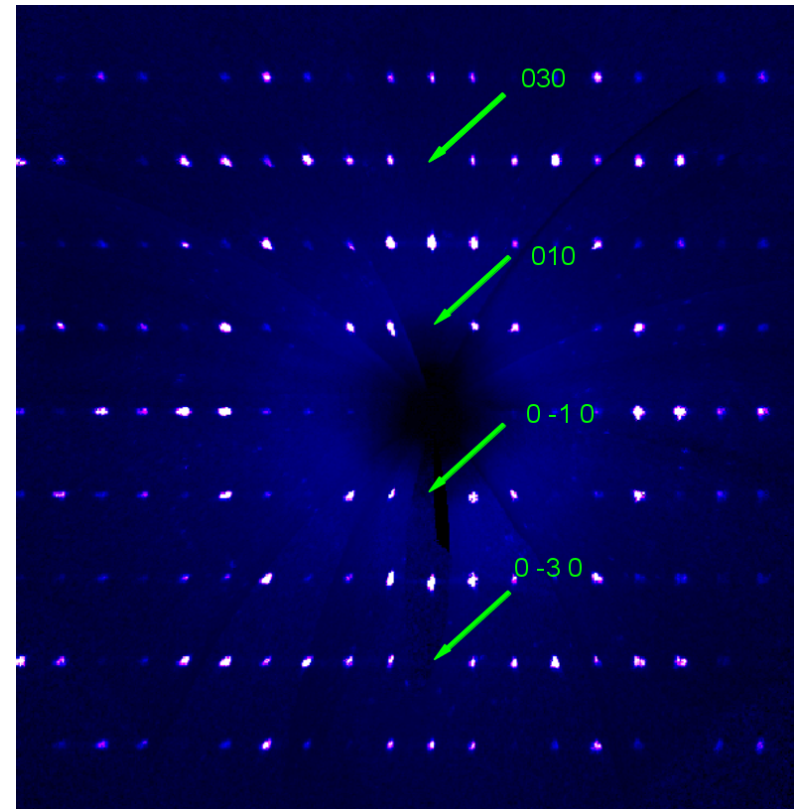

# Orientation of Glide Planes

- When a glide plane is present one can determine the orientation and type of glide plane present from the affected class(es) of reflections
- The 0 index of the affected layer indicates the orientation of the glide's reflection
  - $0kl$ : glide reflects across  $(100)$
  - $h0l$ : glide reflects across  $(010)$
  - $hk0$ : glide reflects across  $(001)$

# Identification of Glide Planes

- The translational component identifies the type of glide plane
- The translational component causes absences in along the affected axes
- $0kl$ :  
 $k = \text{odd} \rightarrow b \text{ glide}; l = \text{odd} \rightarrow c \text{ glide}; k+l = \text{odd} \rightarrow n \text{ glide}$
- $h0l$ :  
 $h = \text{odd} \rightarrow a \text{ glide}; l = \text{odd} \rightarrow c \text{ glide}; h+l = \text{odd} \rightarrow n \text{ glide}$
- $hk0$ :  
 $h = \text{odd} \rightarrow a \text{ glide}; k = \text{odd} \rightarrow b \text{ glide}; h+k = \text{odd} \rightarrow n \text{ glide}$

# Example of c glide ( $h0l$ : $l = \text{odd}$ )

